



Product Information

The Modern and Long-lasting
IRPL® Treatment Solution

E-SWIn ^{usa}



Innovative and safe technology

- Engineered, designed and manufactured in France by E-SWIN, a leading manufacturer of high-tech medical equipments.
- First engineered and designed IRPL® for treating the periorbital area.
- Painless, safe, and reliable.
- Using a patented technology (IRPL® - Intense Regulated Pulsed Light) to generate perfectly calibrated light pulses.
- FDA cleared using the 510k predicate for the treatment of rosacea.
- Compact & Portable, no space restrictions, with a ready-to-go travel case for satellite locations.
- Clinically proven safety and efficacy.

Treatment protocol

- The treatment only takes a few minutes per session.
- Patients are comfortably seated on a treatment chair.
- The protective white metallic eyemask is adjusted on the patient's eyes.
- Then, a specific treatment gel is applied on the cheekbone and the temporal areas.
- In total five flashes are applied on each eye. Four flashes are applied in the lower zone of both - starting on the inner canthus - and one flash in the temporal area to stimulate the upper glands.
- The gel is wiped off and the treatment is complete.
- No specific care after the treatment is required.




Instant & cumulative improvements

The **E•Eye** treatment is performed within a short, efficient protocol of 3-4 sessions. The effect of the treatment appears very rapidly after each session and is cumulative. Usually the results last a minimum of 6 months. In order to improve the achieved results and to avoid returning discomfort, a yearly maintenance

session is recommended to continue to achieve optimal result.

This treatment protocol consists of 3 sessions (and a fourth optional session) as follows:

**Day 1 / Day 15 / Day 45 /
Day 75 (optional)**

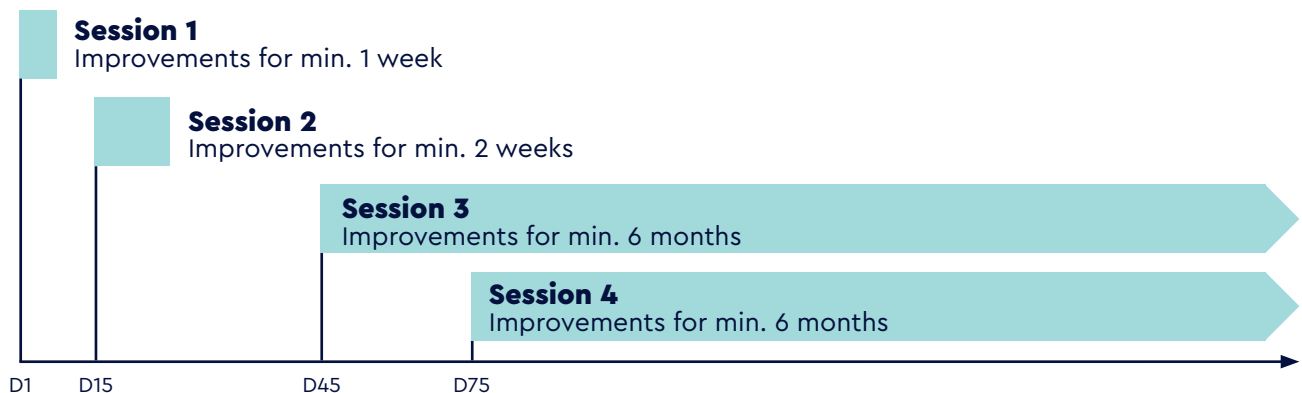

Product specifications



Dimensions	345 × 320 × 440 mm
Weight	11,5 kg
Dimensions Packaging	740 × 460 × 610 mm
Wavelength	> 580-1200 nm
Technology	IRPL® (Intense Regulated Pulsed Light)
Manufacturer warranty	2 years
Device Service	Maintenance-free
Installation type	Desktop unit, mobile use
Registrations	ANVISA, CE, CFDA, FDA, HC, TGA,...



Made by E-SWIN



Medical Certification



Patented Technology

Learn more about our innovations.

Contact your local **E>Eye** specialist

© E-SWIN USA. All rights reserved. Ref. Nr. M029US050B08A rev. B 2023-11-14
E>Eye is FDA cleared under 510(k) K200616 for the following indication for use: E>Eye is a prescription device indicated for the treatment of Rosacea.



Discover the difference.

Smart screening &
Long-lasting Treatment

www.e-swin.us

www.keeptoureyesfit.com